



Pulsair intelliPuff

Non-Contact Tonometer

Pulsair intelliPuff

The Keeler Pulsair intelliPuff Non-Contact Tonometer is indicated for measuring intraocular pressure without contacting the eye to aid in screening and diagnosis of glaucoma.




Key features

- Time on target system – easier to use
- Algorithm to ensure best readings – greater accuracy
- Lighter puff – soft and gentle for your patients
- Space saving flexibility - Pulsair intelliPuff can be wall or desk mounted
- The Pulsair intelliPuff requires no sterile consumables so ongoing costs are managed

Convenient and easy to use

Being hand held, the Pulsair intelliPuff avoids the need to position the patient in a chinrest removing the entire process from the patient journey, making it quick and easy to use while the patient is already seated in an examination chair.

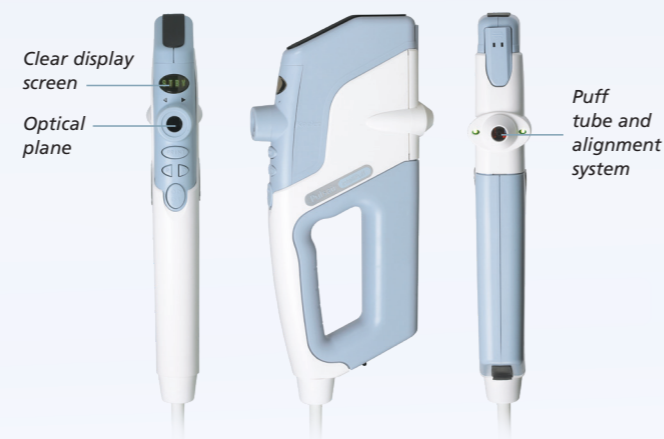

Soft and gentle air puff technology

Specifications

Pulsair intelliPuff

CALIBRATED RANGE	5mmHg to 50mmHg
MEASUREMENT SCALE	mmHg (millimetres of mercury)
DISPLAYED ACCURACY	0.1mmHg
DISPLAY	4 character dot matrix scrolling
FIRING SYSTEM	Automatic and objective
ILLUMINATION SYSTEM	LED infra red
WORKING DISTANCE	20mm
MOUNTING SYSTEM	Table/desktop/wall
UMBILICAL CORD LENGTH	2.0m

Dimensions and Weight

BASE UNIT	260 x 220 x 215mm (HxWxD) Weight 2.465Kg
HANDPIECE	315 x 46 x 150mm (HxWxD) Weight 0.890Kg

Printer

Thermal line printer



Space saving - the Pulsair intelliPuff can be wall mounted

Thermal line printer



Part Numbers

Pulsair intelliPuff Non-Contact Tonometer
2414-P-2001 Pulsair intelliPuff Tonometer

Accessories

2208-L-7008 Printer paper roll
2414-P-7011 intelliPuff wall mounting kit
EP39-70304 Pulsair intelliPuff dust cover
2415-P-7038 Tonometer face shield



Pulsair intelliPuff Overview

The Pulsair intelliPuff is a portable screening device used to measure intraocular pressure which is one of the key biomarkers considered in diagnosing glaucoma.

The Pulsair intelliPuff is an electrically powered handheld type Non-Contact Tonometer. The system provides the user with a position detection display to aid alignment of the tonometer relative to the eye.

Once the system is aligned, the Pulsair intelliPuff will automatically take IOP measurements. The system also allows the user to print IOP measurements taken of both eyes.

What's in the box?

- Pulsair intelliPuff
- Power Supply Unit
- Pulsair intelliPuff Dust Cover
- USB Instruction Manual
- intelliPuff wall mounting kit
- Thermal printer roll
- Tonometer face shield

More Keeler Instruments

With over 100 years of industry experience, Keeler has become a household name among optometrists and ophthalmologists worldwide.

View our full range at

www.keeler.co.uk | www.keelerusa.com

Got a question? - Call us

UK: +44 (0)1753 857177 India: +91 99 3031 1090

USA: +1 610 353 4350 Brazil: +55 11 4302 6053

China: +86-18512119109

Keeler
— A world without vision loss —



A **Halma** company

Manufacturing address



Keeler Ltd. Clewer Hill Road Berkshire SL4 4AA UK

CE 1639